

Rigi-Box®

steri

Safety box

for sterilization and disposal
of medical waste.

(vacuum-steam-vacuum at 134 °C)




H steri 35 = 33 cm / H steri 60 = 52 cm / L x B steri 35 und 60 = 39 cm x 39 cm

Types: **Rigi-Box steri**
35 and 60

material:
polypropylene

The box have a
stable handle in
the lid and two
integrated handle-
hollows at the
bottom part.

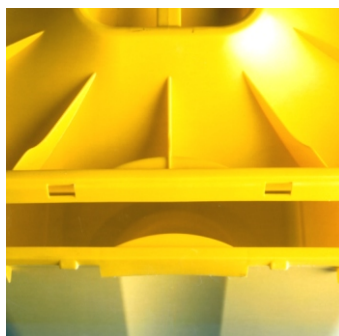

A safe and stable box for
disposal of A- and B-waste
or infectious C-waste, which
will converted into B-waste
through the sterilization in the
sterilization equipment.



The little red lid refines
at about 110 °C and let
the opening for the
sterilization process
automatically through
(vacuum-steam-
vacuum).
The waste can not fall
out because of the
bars which are in the
lid.



The box can be
stacked in empty
state as well as in
full state and is
therefore space-
saving.

Deutsches Patent
Nr. 197 32 979.9 / Nr. 199 01 269.5
UN-Nr.: **Rigi-Box steri 35**
UN/3H2/Y12/S/*/B/RIG-DID-000039
UN/3H2/Y12/S/*/B/RIG-DID-000040
UN-Nr.: **Rigi-Box steri 60**
UN/3H2/Y15/S/*/B/RIG-DID-000043

Medical Disposal Containers

Industriestraße 44 • D-75382 ALTHENGSTETT

Phone +49 (0)7051 1745 • Fax +49 (0)7051 7418

Internet <http://www.rigling.com>

eMail rigling@t-online.de

Rigling

GmbH